

WIRED NURSE CALL PLUG TYPE GUIDE

When you order a wired Alert Mat or device, choose the product code ending that matches your nurse call system. That way, you'll get the right plug.

Don't go by plug type alone. Some plugs look identical but are wired differently inside to suit specific nurse call systems. If you order the wrong code ending, the product may not work properly.


Use this guide to match your nurse call system to the correct product code ending.

Nurse Call System	Product Code Ending	Plug Type
Aid Call ATX3000 (old)	-W	White 'reverse BT type' plug
Aid Call ATX3000+	-C6TSPPro	RJ12 clear plug, 6P6C
Aid Call ATX4000	-C6	RJ12 clear plug, 6P6C
Aid Call ATX5000	-C6TSPPro	RJ12 clear plug, 6P6C
Aid Call Touchsafe Pro	-C6TSPPro	RJ12 clear plug, 6P6C
ARM (Alarm Radio Monitoring)	-D	Double ring stereo jack plug
Arquilla	-C6TSPPro	RJ12 clear plug, 6P6C
Chubb	-S	Single ring mono jack plug
Chubb	-C6-Chubb	RJ12 clear plug, 6P6C
CLB Nurse Call	-CLB (wireless only)	3.5mm 4 pole jack plug
Courtney Thorne	-C	Single ring mono jack plug
C-TEC Nursecall	-Q	Double ring stereo jack plug
Eclipse	-D	Double ring stereo jack plug
Gemini	-SASN2 (wireless only)	Double ring stereo jack plug
Intercall	-S	Single ring mono jack plug
LiveLink	-LL	RJ9 clear plug, 4P4C
Medicare (floor mats)	-MED01	Clear 'telephone' plug
Medicare (bed/chair pads)	-MED22	Clear 'telephone' plug
One Call Systems	-S	Single ring mono jack plug
Quantec	-Q	Double ring stereo jack plug
SAS NetworkII	-SASN2 (wireless only)	Double ring stereo jack plug
Stanley Blick	-SASN2 (wireless only)	Double ring stereo jack plug
Wandsworth	-S	Single ring mono jack plug

Plug Type Identification



Single ring mono jack plug

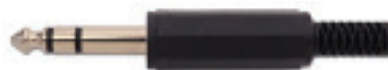

Double ring stereo jack plug



RJ9 Clear Plug, 4P4C



RJ12 Clear Plug, 6P6C



White 'reverse BT type' plug