

# MARSDEN

## M-610

### Marsden Portable Wheelchair Weigh Beams



Class III Approved scale for weighing patients in wheelchairs

Portable with handles; carry cases available

Set up and weigh your patient in seconds

Accurate to 100g

[www.marsden-weighing.co.uk](http://www.marsden-weighing.co.uk)



## Marsden M-610: An accurate weight for patients in wheelchairs - in seconds

Marsden's M-610 means weighing patients in wheelchairs is easier than ever - because you can weigh them accurately, **without them even needing to leave their wheelchair**. The M-610 is **accurate**, **portable** and **setting up is easy**.

### Quick and easy set-up

At Marsden, we believe that a wheelchair should not be a barrier to getting an accurate weight reading of a patient.

Therefore, the M-610 is accurate, easy to use and has been made completely portable. The two beams have carry handles, and once in position, fold down the ramps, and then the beams simply need to be connected to the **premium DP-3810 indicator**.

Once the scale is switched on, simply make sure the beams are the correct distance apart, and weigh your patient.


### Calculate patient BMI automatically - and then record it


If you already know the height of the patient, you can enter this using the **full numeric keypad**. The scale will then **calculate the patient's BMI**.

This means you'll have an accurate BMI reading - which then can be printed, along with the patient's weight, using the **optional TP-2100 Printer**.

Additionally, choose **Bluetooth or Wifi connectivity** and send this information **directly to a patient database**.

[www.marsden-weighing.co.uk](http://www.marsden-weighing.co.uk)

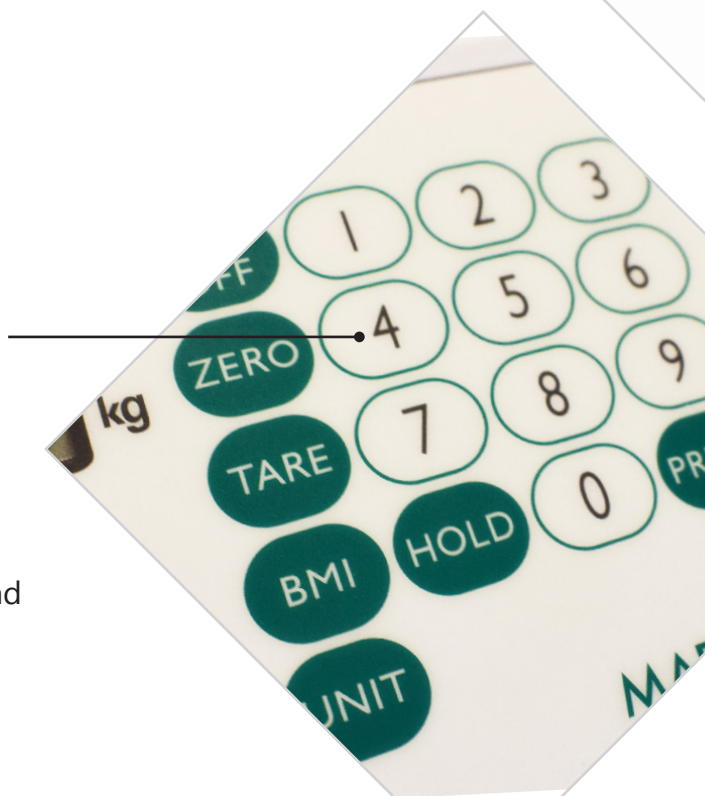

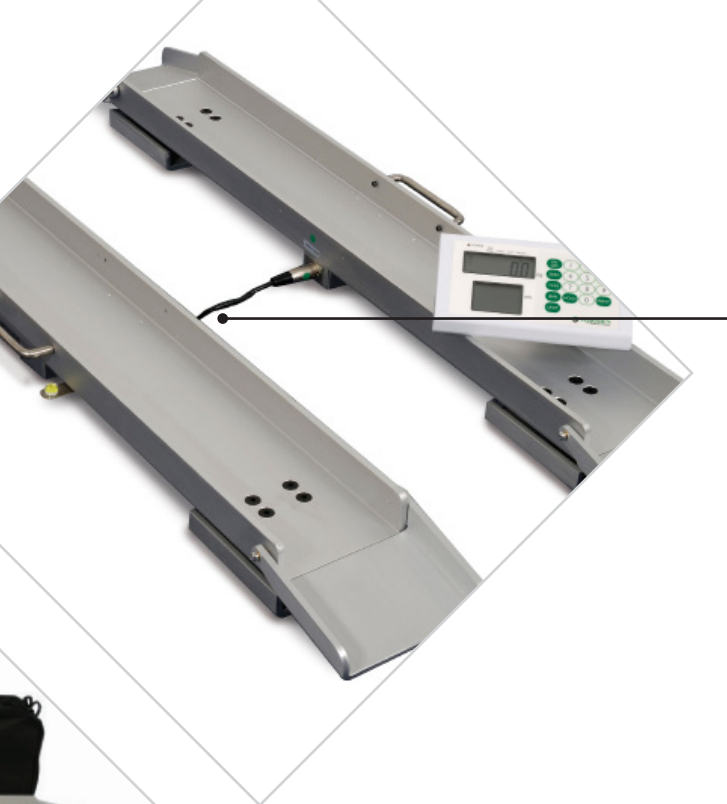


## Less wires, less mess

The M-610 is not only easy to transport - there are **carry cases available** too - it also has an **internal rechargeable battery**.

An adaptor is supplied, and fully charged, the battery will provide a minimum of **40 hours' continuous use**.

## Key features

- Class III Approved
- Rechargeable battery (adaptor provided)
- Full numeric keypad
- Three displays - weight, height, BMI
- BMI calculation
- Hold/Tare/Zero
- Low profile ramps
- Lightweight, portable with carry handles
- 40 hours' continuous use from full charge
- Optional printer
- Optional carry cases
- Optional Bluetooth, Wifi or Micro USB

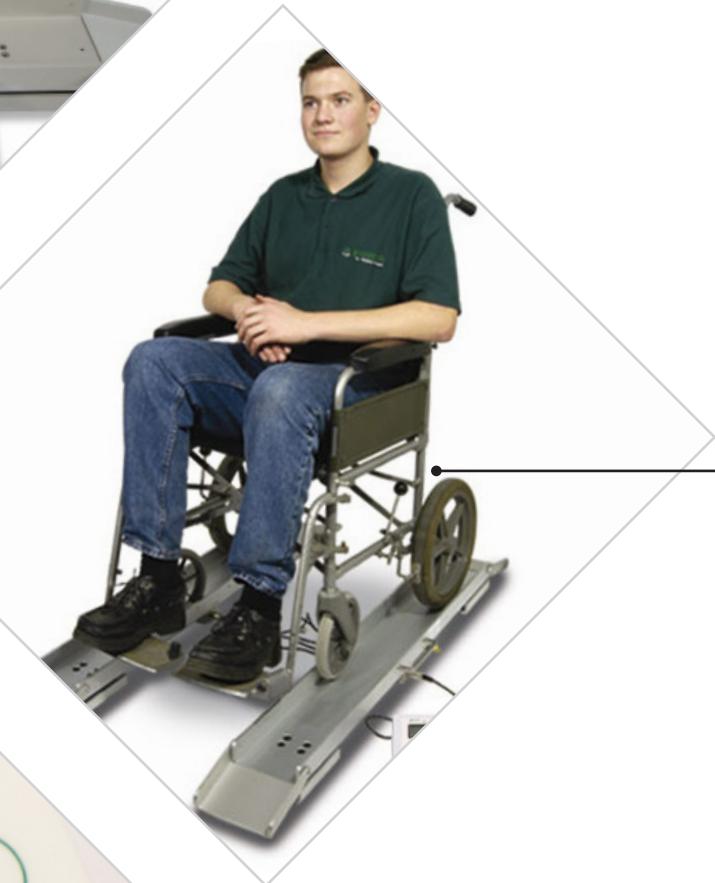

## A capacity suitable for most patients

A **300kg capacity** means the M-610 is suitable for almost any patient. It'll weigh to a **medically-approved accuracy of 100g**.

If you wish to exclude the weight of the wheelchair from your readings on the scale, you can do this using the **Tare function**.


## Specification

- Class: III
- Capacity: 300kg
- Graduations: 100g
- Weight of Beams: 8.5kg
- Power: Mains adaptor/rechargeable battery
- Beam dimensions: 1200mm x 150mm x 95mm

[www.marsden-weighing.co.uk](http://www.marsden-weighing.co.uk)



# Marsden M-610: Key features at-a-glance

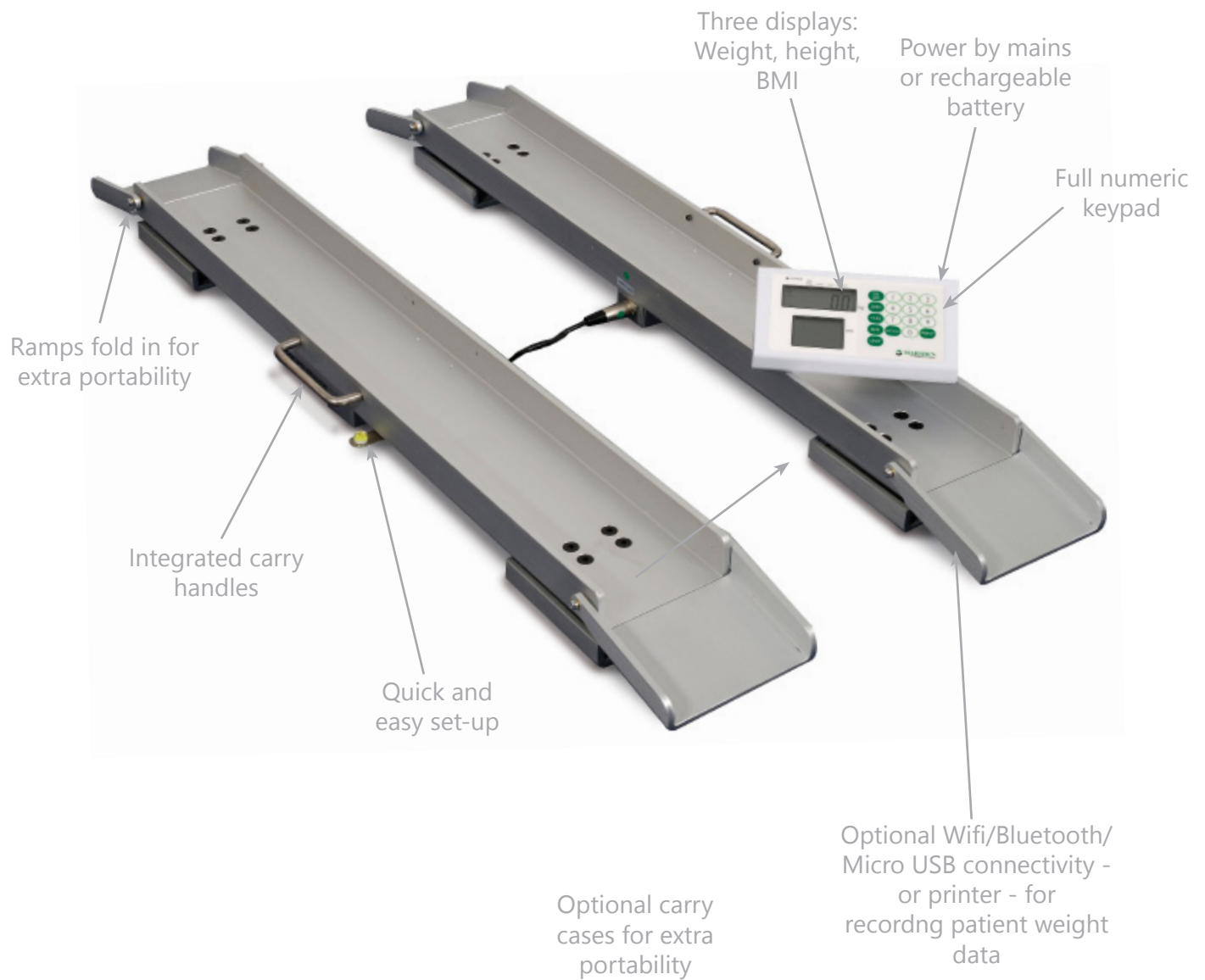

# MARSDEN

Tel: 01709 364296 / 0800 169 2775

Fax: 01709 364293

E-mail: [sales@marsdengroup.co.uk](mailto:sales@marsdengroup.co.uk)

[www.marsden-weighing.co.uk](http://www.marsden-weighing.co.uk)

 MarsdenWeighing

 @MarsdenWeighing

 [www.youtube.com/MarsdenWeighing](http://www.youtube.com/MarsdenWeighing)

 <http://pinterest.com/marsdenweighing>

