

ECTO(VISE

Face Visor



- Anti-fog
- Comfy foam with elastic headband
- Visor suitable for use with prescription lenses and protective masks
- Can be used for prolonged period – does not create distortion or fatigue
- Part of recommended PPE for protection from viruses
- Medically approved
- Manufactured in the UK in quality controlled, clean production areas
- Individually packed

CE




 **PolarSeal**

polarseal.net

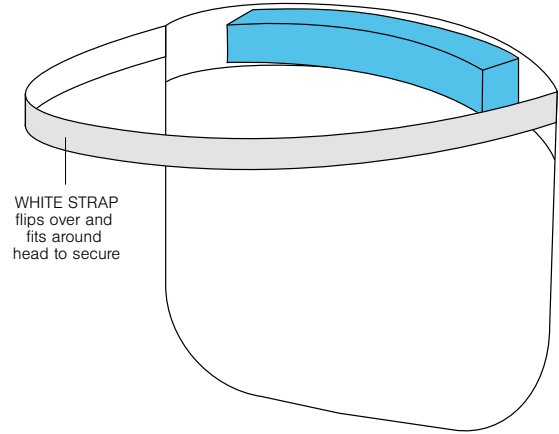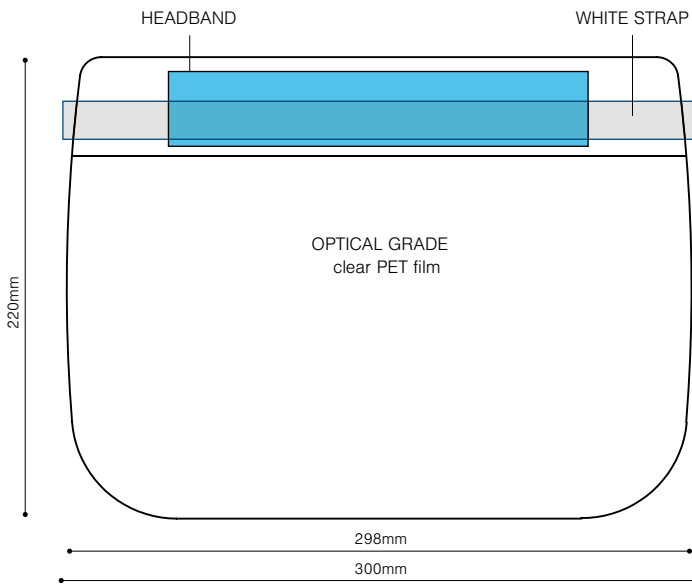

1 Wash hands thoroughly for 20 seconds

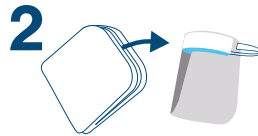

2 Carefully remove ECTO(VISE from packaging



3 Place over face and position strap around back of head



4 Ensure foam is on forehead and face is fully screened

Specification	Detail
PRODUCT CODE	ECTO(VISE
PACK QUANTITY	25 or 50
PRODUCT COLOUR	PERFECTLY CLEAR (NO TINT) – SURGICAL STANDARD
PRODUCT TYPE	FACE SHIELD
OVERALL SIZE	300mm (W) X 220mm (H)
WIDTH OF VISOR ONLY	298mm
HEIGHT OF VISOR ONLY	220mm
FIELD OF USE	LIQUIDS AND DROPLETS BARRIER
APPLICATION	MEDICAL, GENERAL
STANDARDS	EN 166:2001, CE MARK

Material	Description	Size
VISOR FILM	OPTICAL GRADE PET	0.25mm THICK
FOAM HEADBAND	MEDICAL GRADE SOFT GRADE PU	30mm THICK
FOAM HEADBAND ADHESIVE	PSA	N/A
ELASTIC HEADBAND	MEDICAL GRADE WOVEN ELASTIC	19mm THICK

Additional benefits and details

- Medical grade materials
- Single use
- Can be used for prolonged periods as it does not create distortion or fatigue
- Visor suitable for use with prescription lenses and surgical masks
- Always check before using